

Open Base Power XL Table with Adjustable Backrest

MODEL 80299

- Features easy-clean, 30" wide top with adjustable pneumatic backrest which is adjustable from either side and raises 80 degrees from flat
- 600 lbs. (272 kg) load capacity under normal use
- 5-year limited warranty on frame and electric motor drive
- Smooth, electric, raise and lower height adjustment
- Pneumatic backrest cylinder is mounted internally
- Full width backrest hinge adds rigidity for greater support
- Heavy duty, steel frame is fabricated in the USA by Clinton Industries
- Dependable scissors-lift features an H-leg design for extra stability
- Rounded corner 3/4" plywood top with 2" of ultra firm padding covered by heavy weight upholstery
- Harden steel bolts and self lubricating pivot bearings for years of dependable service
- Low height may help achieve ADA tax credit

CLINTON OPEN BASE & SHROUDED
POWERXL
Extra width, Extra capacity, Extra Affordable




80299 SPECIFICATIONS

| | |
|---------|--------------------------|
| Overall | 72" L x 30" W x 18-35" H |
|---------|--------------------------|

| | |
|----------|---------------|
| Backrest | 38" L x 30" W |
|----------|---------------|

California residents please be advised of the following, pursuant to Proposition 65:

1
 **WARNING:** This product can expose you to chemicals, including chromium, which is known to the State of California to cause cancer or reproductive toxicity. For more information go to www.P65Warnings.ca.gov.



Upholstery Colors



Slate Blue



Wedgewood



Royal Blue



Desert Tan



Burgundy



Country Mist



Warm Gray



Gunmetal



Black

OPTIONS

Safety Straps Set

- 2" x 92" strap is secured with quick hook and loop Velcro® fasteners
- 2 safety straps per set (92" x 2")

080



Safety Rails (pair)

- 3 height positions
- Plunger height adjustment

082



Hand-Held Controller

- Alternative to foot control
- Stores easily under top when not in use

084



Locking Casters

- Heavy duty, rubber wheel, easy-roll casters

087



030 Paper Dispenser

- Accommodates paper rolls from 18" (45.72cm) to 21" (53.34cm) wide


030



Paper Cutter

- Clear plastic with elastic straps and buckle

040

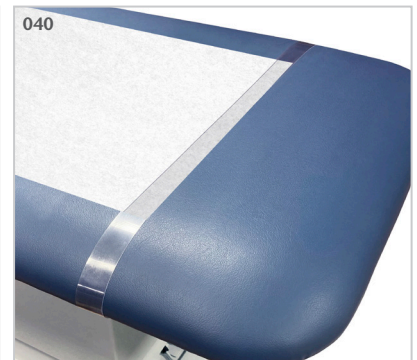

Paper Dispenser/Cutter Combo

- Paper dispenser and paper cutter ordered as one option

034