

## Premium Guard Gel/Foam Overlay

14893 Gel Overlay, 76" (L) x 34" (W) x 3.5" (H), 300 Wt. Cap., 1/cs  
 HCPCS code: E0185  
 14901 Bariatric Gel Overlay, 76" (L) x 42" (W) x 3.5" (H), 500 Wt. Cap., 1/cs  
 HCPCS code: E0185

### FEATURES:


- New and improved Gel Overlay has an upgraded **300 lb. weight capacity** and optimizes pressure redistribution.
- New upgrade: 10 gel bladders cover entire patient contact area.
- Full length "Head to Toe" Gel Bladders minimizes "bottoming out" and "pitting".
- New convoluted, egg-crate top layer vastly improves patient comfort.
- 4 Elastic corner straps secures overlay to mattress.
- Urethane coated nylon cover is, water resistant and vapor permeable.
- Vinyl base is stain resistant and waterproof.
- Bladders contain an aqueous non-toxic gel.
- Designated head and foot sections.
- 1 Year Limited Warranty on cover and mattress.

### SURFACE COMFORT TOP LAYER

### GEL LAYER

### GEL OVERLAY

### PREMIUM GUARD OVERLAY WITH SOFT NYLON COVER!

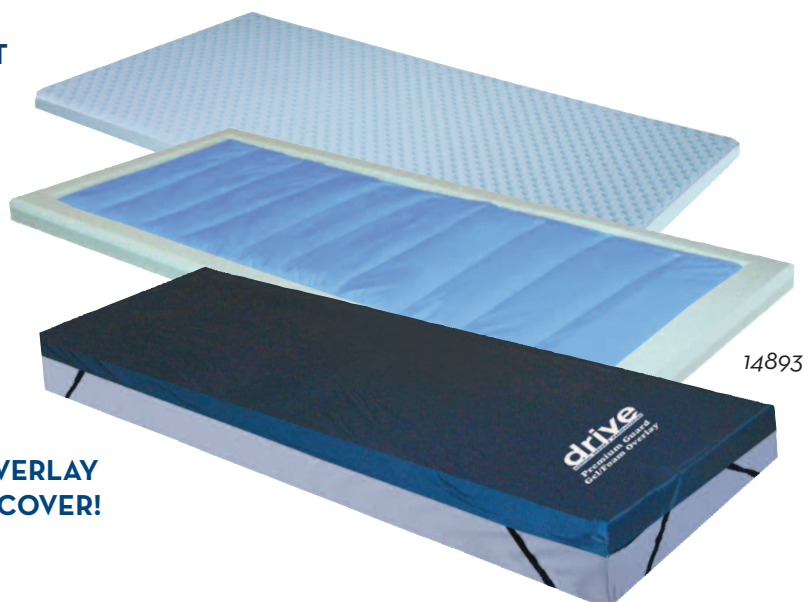

14904

### Bed Side Rail Safety Pads

14904 70" (W) x 13" (D) x 1" (H), 1pr/cs

### FEATURES:

- Prevents semi-conscious, comatose and overdosed patients from injuring their limbs against the bed side rails.
- Manufactured with foam core.
- Covered in vinyl with velcro straps to securely fasten to bed side rails.
- Folds in three sections for convenience of use and storage.

### SPECIFICATIONS

Item #	Length	Width	Height	Product Weight	Carton Shipping Weight	Weight Limit
14893	76"	34"	3.5"	47.8 lbs.	50.8 lbs.	300 lbs.
14901	76"	42"	3.5"	50.5 lbs.	53.5 lbs.	500 lbs.
14904	70"	13"	1"	2.3 lbs.	3.3 lbs.	NA

\*Important Notice: Item numbers that display an asterisk (\*) reflect an item that will be discontinued or package change.

pressure prevention

To order call toll free: 877.224.0946

drive  
WWW.DRIVEMEDICAL.COM