

OUR MOST POPULAR MODEL JUST BECAME...

EVEN BETTER!

DELTA™

Ultra Light 1000, Semi-Electric Bed



**New & Improved
Game Changer**

BED, SIDE RAILS AND MATTRESS SAFETY-TESTED BY SGS, A NATIONALLY RECOGNIZED TESTING FACILITY, SATISFY ALL APPLICABLE FDA ENTRAPMENT GUIDELINES.

15030

beds

To order call toll free: 877.224.0946

drive™




Truly universal with most head and footboards by other manufacturers. This makes managing inventory much easier!



Improved rail placement. Labels and color-coded springs assure accurate installation and patient safety.



Head-board taller than foot-board. Better fits home decor!

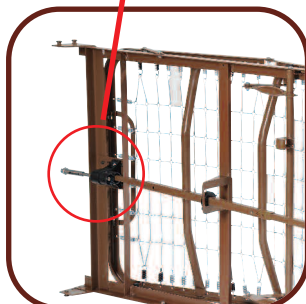

Transition box is on the foot section NOT on the head and footboards. That's key to its interchangeability!



High-low shaft mounted on bed. Once installed, it can't be lost or forgotten on future installations!

Important Notice: Item numbers that display an asterisk () reflect an item that will be discontinued or package change.

DELTA™ Ultra Light 1000, Semi-Electric Bed

Available
1 Pack

15030	Ultra Light Semi-Electric Bed
15030-PKG-11	Ultra Light Semi-Electric Bed with Full Length Side Rails (item #15001ABV) and 80" Quick 'N Easy Comfort Mattress (item #15076)
15030-PKG-12	Ultra Light Semi-Electric Bed with Half Rails (item #15201BV) and 80" Quick 'N Easy Comfort Mattress (item #15076)
15030BV-FR	Ultra Light Semi-Electric Bed with Full Length Side Rails (item #15001ABV),
15030BV-HR	Ultra Light Semi-Electric Bed with Half Rails (item #15201BV),
15030BV-PKG	Ultra Light Semi-Electric Bed with Full Length Side Rails (item #15001A/15001ABV) and 80" Inner Spring Mattress (item #15006)
15030BV-PKG-T	Ultra Light Semi-Electric Bed with Full Length Side Rails (item #15001ABV) and 80" Therapeutic Support Mattress (item #15019)
15030BV-PKG-1	Ultra Light Semi-Electric Bed with Half Rails (item #15201BV) and 80" Inner Spring Mattress (item #15006)
15030BV-PKG-1T	Ultra Light Semi-Electric Bed with Half Rails (item #15201BV) and 80" Therapeutic Support Mattress (item #15019)
15030BV-PKG-2	Ultra Light Semi-Electric Bed with Full Length Side Rails (item #15001ABV) and 80" Foam Mattress (item #15007)
15096-SEMI	Ultra Light Semi-Electric Bed Reclamation Kit, 1/ea (Includes sleep surface & motor only, no head and foot boards)

FEATURES:

- The Delta 1000 is truly universal because the headboards and footboards are completely interchangeable with Drive's and most other manufacturers bed ends, old and new.
- The transition box mounted on the foot section is key to interchangeability.
- The transition box allows you to change the rotation of the high-low shaft.
- Headboard is taller than footboard to better fit home decor.
- The wood grain panels are break and scratch resistant.
- Foot section weighs 50% less than conventional foot sections to make deliveries safe and easy.
- Once installed, the high/low shaft will stay with the bed to eliminate loss or forgetting on future installations.
- Labels and color-coded springs assure accurate installation of side rails for patient safety.
- Unique motor is completely self-contained to reduce weight and noise.
- The motor assembly can be installed and removed with patient in bed.
- The motor is UL approved.
- In the event of a power failure, a 9 volt battery, located in the motor, can lower the head and foot sections nine times. No crank needed.
- New and improved hand pendant has large, easy to use controls.
- Bed ships in two cartons.
- Channel frame construction provides superior strength and reduced weight.
- Zinc coated spring deck.
- Extension kit available model # (15030EXTKIT) extends length from 80" to 84".
- **Warranty: (Original purchaser)**
5 years on motor and hand control, Limited Lifetime on welds and frame,
1 year on all other parts and components.

SPECIFICATIONS

Item #	Bed Height (min) w/wheels	Bed Height (max) wo/wheels	Overall Length	Overall Width	Product Weight	Carton Shipping Weight	Weight Limit
15030	12.75"	24"	88"	36"	155 lbs.	170 lbs.	450 lbs.

Important Notice: Item numbers that display an asterisk () reflect an item that will be discontinued or package change.