

INVACARE®
APM™
MATTRESS REPLACEMENT SYSTEM



The **Invacare® APM™ Mattress Replacement System** is an easy to use alternating pressure mattress that aids in the prevention and treatment of Stage I to Stage IV pressure ulcers. This economical system does not compromise quality with price and is reliable with low maintenance.

Features

State of the art compressor – With solenoid valve system allows for easy set up and effective therapy cycles

Therapy time – Provides patient comfort and relief with 10 minute cycles

Constant pressure monitoring – Control unit constantly monitors the air pressure in the mattress to maintain the appropriate pressure to prevent capillary occlusion (generally below 32mm/Hg)

Pressure adjustment – Pressure can be adjusted by 25%, 50% or 75% to provide comfort for frail patients

Static mode – Alternating function ceases and pressure is redistributed across all cells within the mattress

Max firm mode – Assists during patient ingress/egress and initial set up

Audio and visual alarms – Warns patient and caregiver during power outage or low pressure

Power outage safety – Mattress remains inflated during power outage as long as air hose is connected to the compressor

HCPCS Code: E0277



- 1. Digital display with settings based on patient weight
- 2. Audio and visual alarms
- 3. Quilted cover with 600 denier nylon bottom




Invacare Corporation
www.invacare.com

USA
One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada
570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

©2012
Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.
Form No. 10-166
Rev. 03/12
111231

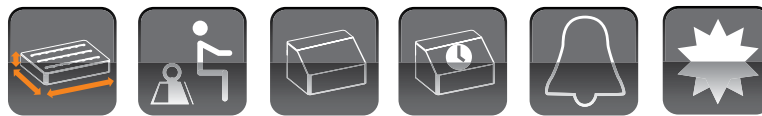

Model #	Dimensions	Weight Capacity	Compressor	Alternating Time Cycle	Alarms	Limited Warranty
CG9900	80" L 36" W 6" H	350 lb.	10 LPM	10 min.	Audio and Visual	1 Year (Mattress) 1 Year (Power Unit)

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.

