


Cap-Shure® Cotton Tipped Applicator

Polystyrene handle

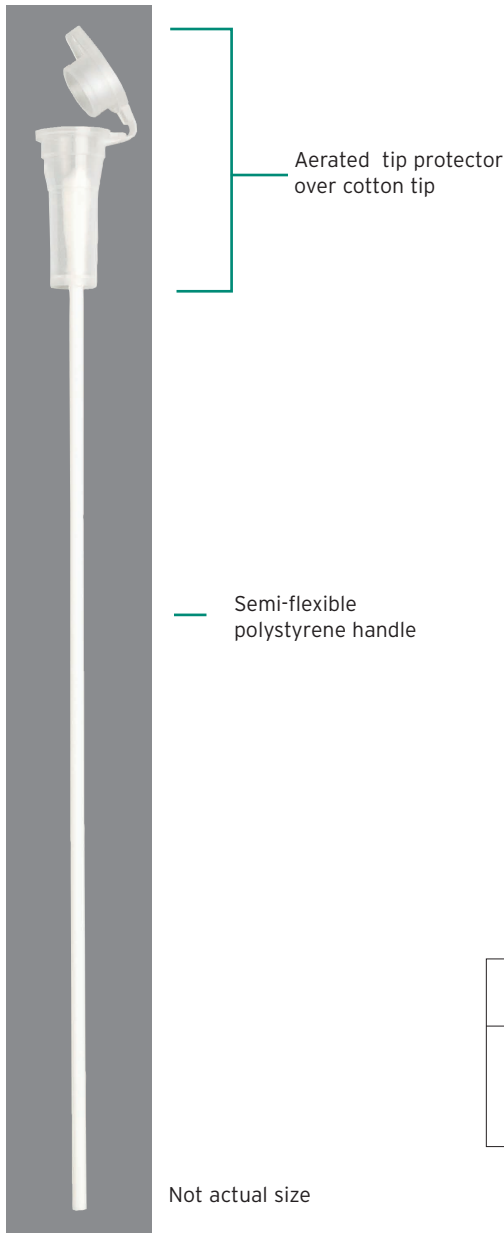

Description

This absorbent cotton tipped applicator is produced with High Quality USP grade cotton which is carded and coiled onsite to maintain consistency. It is bonded securely to semi-flexible polystyrene handle by an aqueous based adhesive. Aerated tip protector designed for convenient, reliable specimen, DNA or evidence collection. Comes as a unit with a reclosable, ventilated cap for air drying and transporting samples.

Dimensions

7" Overall Length (177.8 mm)

Tip

■ 0.203" Dia. (5.156 mm)

■ 0.687" Length (17.45 mm)

Handle

■ 0.100" Dia. (2.54 mm)

■ 6.875" Length (174.625 mm)

Standards

■ CE marked for medical use

■ ISO 9001:2008

■ ISO 13485:2003

■ FDA GMP's

■ Certificate of sterility and/or compliance available with each shipment upon request

Qualities

■ Medical grade quality

■ Tip suspended away from sides

■ Can lay flat to dry

Applications

Diagnostic / Medical / Forensics

■ Specimen collection

■ DNA collection

■ Evidence collection

Ref #	Packaging
25-807 1PC EC <i>Sterile</i>	10/50/1 (10 boxes per case; 50 individually wrapped packages per box; 1 per package)



Quality Medical Products Since 1919

Puritan Medical Products Company LLC
Guilford, Maine 04443-0149 USA

Tel 800-321-2313 (US and Canada Only)
207-876-3311

Fax 800-323-4153
207-876-3130

Web www.puritanmedproducts.com

Email sales@puritanmedproducts.com

1-800-321-2313
puritanmedproducts.com